
MOTEUR C MOTOR C MOTOR C MOTOR

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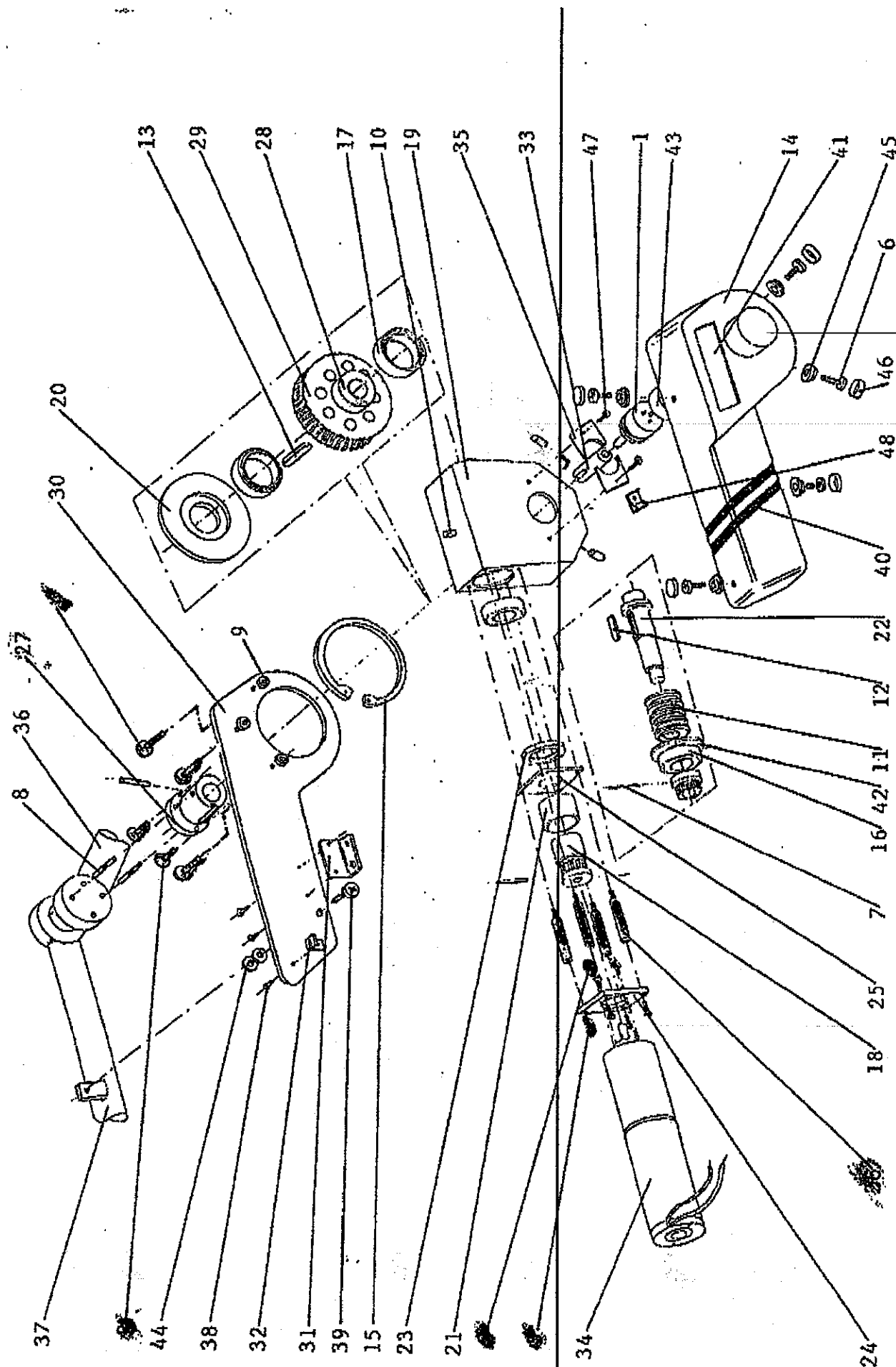COGEMO reserves the right to effect technical modifications.

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
Fax No: (815) 654-8821

Subject: Set for 4071 weight
flex + Ext

STOP POINTS (PULL)
Flex 20 Kilograms
Extend 10 Kilograms
COW - together
weight

KEYPAD FOR 4071 TAW
4091
9461-000-5719
SMITH & NEPHEW ROLYAN INC.
272
FOR DETROIT - GARY

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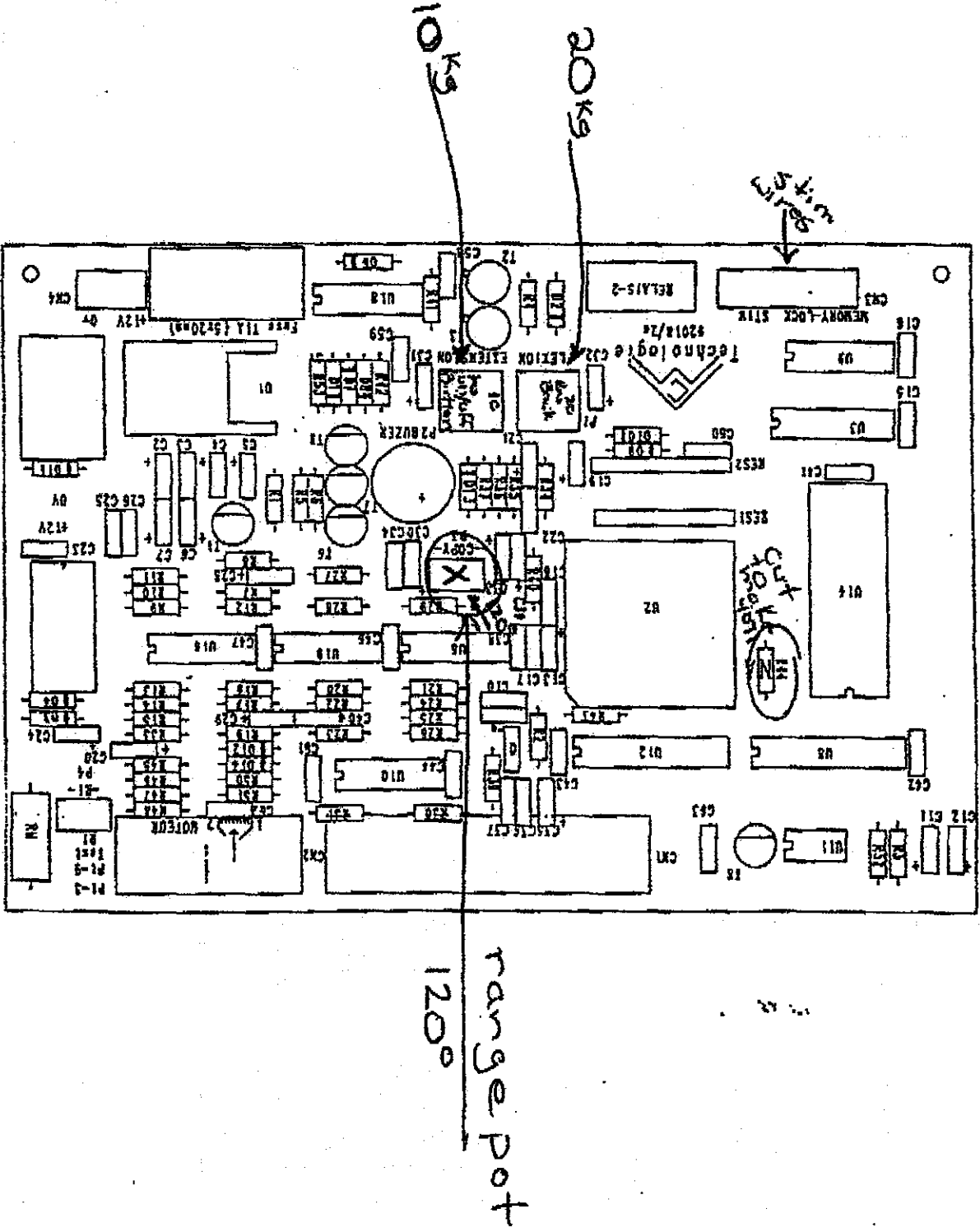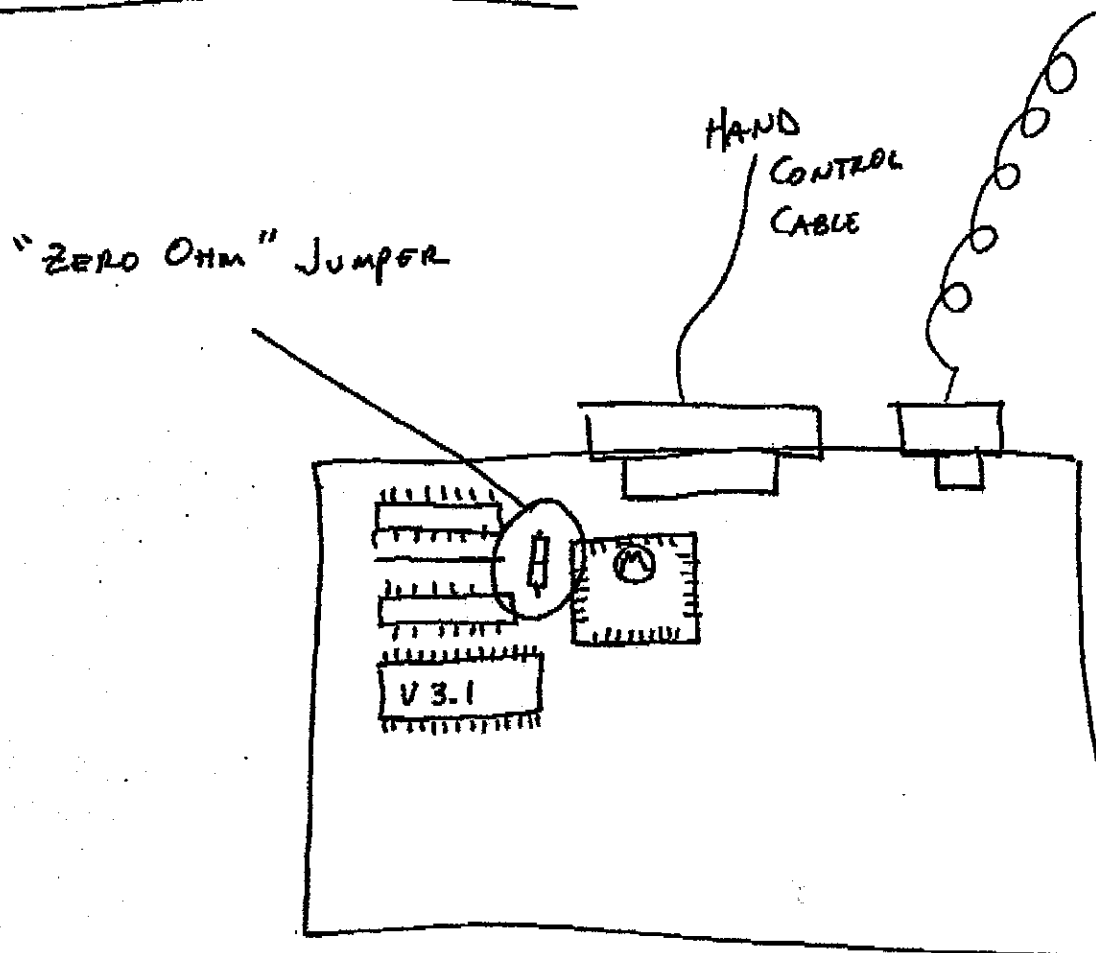
IT IS GREY WITH 1 BLACK STRIPE.

UNDER THE JUMPER ARE THE LETTERS NA.

USE AN OHM METER TO TEST. SHOULD BE LESS THAN 1 OHM. CUT THE JUMPER LEADS AND REMOVE. BOARD IS NOW 4071.

John

VENDOR _____	SHIP TO _____
_____	_____
_____	_____
ATTN: _____	ATTN: _____

Originator: _____ Date: _____

☐ Sale/Resale Sale # _____ Date Order Received _____

Hospital Name or # _____ P.O. # _____

<input type="checkbox"/> Immediate Rental Use <input type="checkbox"/> Departmental <input type="checkbox"/> Credit Memo <input type="checkbox"/> Exchange <input type="checkbox"/> Repair <input type="checkbox"/> Inventory <input type="checkbox"/> Other (See Comments)	<input type="checkbox"/> Mode of Shipment <input type="checkbox"/> UHS Delivery <input type="checkbox"/> Customer P/U <input type="checkbox"/> U.S. Mail <input type="checkbox"/> Special Delivery - Mail <input type="checkbox"/> Ups <input type="checkbox"/> Truck <input type="checkbox"/> Bus <input type="checkbox"/> Air Freight	R.G.D.M # _____ Date _____ Purchaser's Initials _____ <input type="checkbox"/> Ship from Minneapolis <input type="checkbox"/> Drop Ship
-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	-----------------------------------------------------------------------------------------------------------------------------------------------------

COMMENTS: _____

LINE NO.	QTY. ORDER	QTY. SHIP	UNIT	B/C	P/C	VENDOR PART NO.	DESCRIPTION	UHS #	ACCT. CODE	UNIT PRICE	AMOUNT
1						70920432	Bearing ø 9x24x7 mm			6.40	
2						70924075	Axle of Screw with nut and			58.50	
3						70924076	Adjusting washer			3.98	
4						70924078	Press plate of bearing			18.28	
5		4	pc			70920442	Spacer			2.94	
6						70924094	Reductor Spacer			2.38	
7											
8											
9											
10											
11						94670011029	upper case			279.96	
12						94626000852	Nylon pin			.80	
13						946250000671	Covering handle			46.88	
14						946260002701	Washer 4x16mm (for handle)			.80	
15											
16											

Serial Nos. _____	Mode of Shipment (Incoming) _____	Freight (Incoming) Code _____	TOTAL AMOUNT ▶
Reference Nos. _____			Authorizing Signatures:
Date Shipped _____	Mode of Shipment (Outgoing) _____	No. of Pkgs. _____	1. _____
Charges _____			2. _____
Freight (Outgoing) Code _____			3. _____
RECEIVED BY: _____			

DELTA BOOK A
U.S.A.
NO-821-5700

TO: UNIVERSAL HOSP SERVICES
1250 NORTHLAND PLAZA
3800 N 80TH ST
BLOOMINGTON MN 55431

THIS IS AN EXPEDITED ORDER 10107590

SMITH + NEPHEN RICHARDS INC.
4287 DELP STREET, BOX A
MEMPHIS TN 38118 U.S.A.
1-800-821-5700

SHIP TO: UNIVERSAL HOSP SERVICES
5271 ZENITH PKWY
ROCKFORD IL 61111

** PKG. ID # A30515641 **


A305156415111000

11737

SALES REPRESENTATIVE		CUSTOMER ORDER NUMBER		ORDER DATE	
NAME	PHONE	NAME	PHONE	DATE	TIME

Q. LOCATION	CATALOG NO	U/M	QUANTITY ORDERED	QUANTITY SHIPPED		QUANTITY B.O.	DESCRIPTION
				PACKAGES	UNITS		
REFRIG 004550	EA	1	1			0	THIS IS REPAIR ORDER - 000315. GENERAL SURGERY CUSTOMER REPAIR 210058 S/N 505982 S/N 505982 - NO CHARGE UNDER WARRANTY KNURL INSIDE AND REGUE NYLON COUPLING RESET COPY POTENTIOMETER FILE BURR OFF REDUCTOR BOX AT THE POTENTIOMETER MOUNTING HOLE. SMITH + NEPHEN RICHARDS INC.

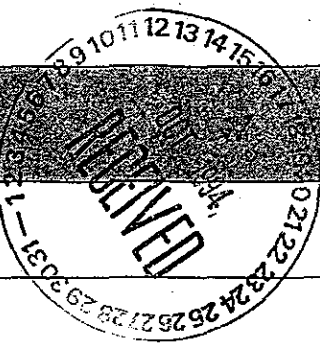Smith + Nephew

4287 DELP STREET, BOX A
MEMPHIS TN 38118 U.S.A.
1-800-821-5700

10-6-94
TAW 0117
SCOTTS BECAUSE:
ANGLES DON'T AGREE
WITH DISPLAY AT
FALL EXTENSION.
JIS

ANY ITEM ORDERED NOT APPEARING ON THIS SHIPPING ORDER
IS BEING SHIPPED FROM ANOTHER LOCATION OR IS ON BACKORDER

1	ZONE 2	ZONE 3	ZONE 4	ZONE 5	ZONE 6	ZONE 7	PACKED BY	TOTAL WEIGHT
G ORDER NO.	ZONE 8	ZONE 9	BOX NO.	BOX NO.	BOX NO.	DATE SHIPPED	POSTAGE/HANDLING	
315641								

I	INTRODUCING THE KINETEC 4071	2
II	CONNECTING THE KINETEC 4071	2
III	HOW TO USE THE HANDCONTROL	4
IV	USING THE KINETEC 4071	10
VI	INSTRUCTIONS IN CASE OF A BREAKDOWN	16

PAGE

- Set up of the pads on the KINETEC 4071
- OPTION : Plastic leg support
- TRANSPORT position
- WORKING position
- Setting up for an adult
- Setting up for a child
- Positionning of the foot plate
- Protection of the opposite foot
- OPTION : Clamp for the foot-end of the bed
- OPTION : Fixation straps (Home use)
- Guarantee limits

- Description : - keyboard
- display
- Locking the handcontrol
- START / STOP function
- Adjusting the EXTENSION limit
- Adjusting the FLEXION limit
- Adjusting the SPEED
- Adjusting the PAUSE
- MANUAL use
- BY-PASSING use
- Timer

- Unpacking the KINETEC 4071
- Safety characteristics
- Starting up the KINETEC 4071

- Definition
- Description


TABLE OF CONTENTS

IMPORTANT
 THE KINETEC 4071 WILL
 NOT RUN IN THE TRANSPORT
 POSITION.
 (SEE PAGE 12)

USER'S
 MANUAL
 4071

Handwritten:
 11/02/92
 J. Legendre

KINETEC®

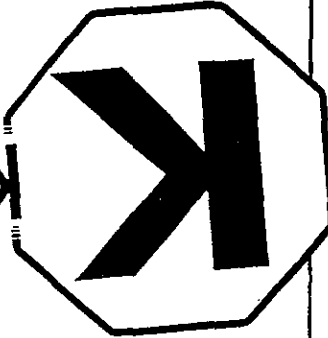Handwritten:
 1800
 821 5700
 DAVE

I INTRODUCING THE KINETEC 4071

Definition

The KINETEC 4071 is a continuous passive motion machine for FLEXION-EXTENSION mobilisation of the hip, knee and ankle joints.

Description

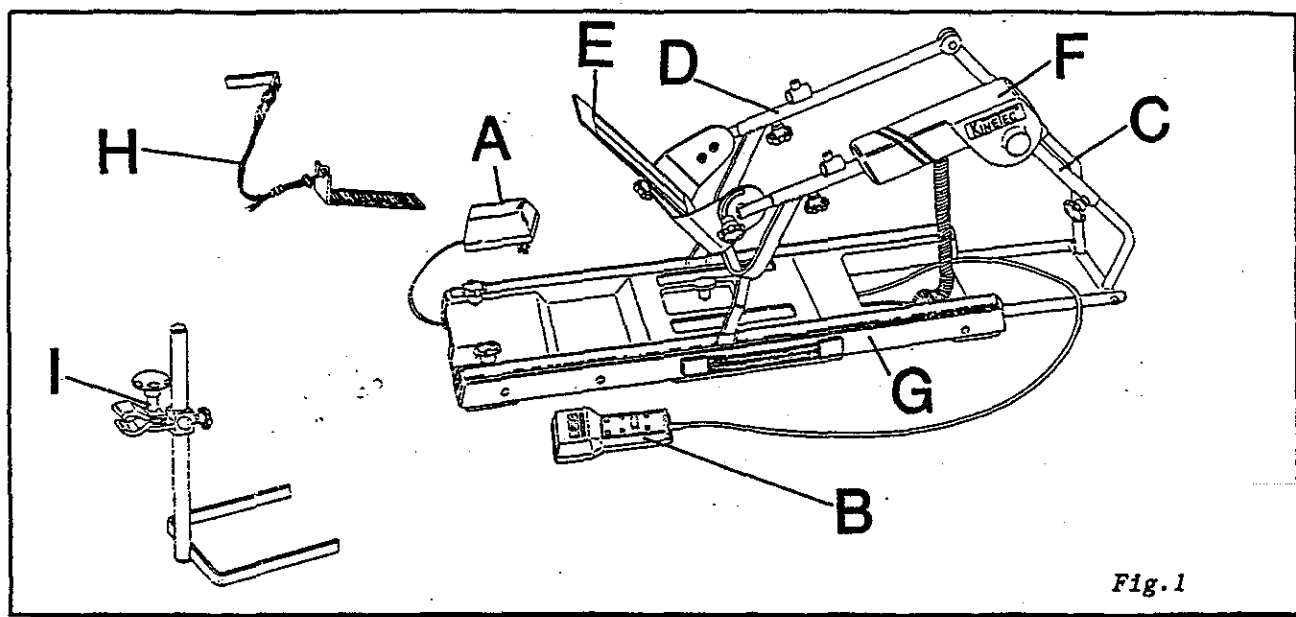
Fig.1


The KINETEC 4071 has the following parts :

- A - Smith + Nephew Richards Power Supply model 7875
(For replacement power supply contact Smith + Nephew Richards
at 1-800-238-7538)
- B - The control box (keyboard and display).
- C - The thigh support, adjustable for length.
- D - The lower leg support, adjustable for length.
- E - The foot plate.
- F - The motor.
- G - The base plate.
It consists principally of the guidance rails and the control
circuitry
- H - Home bed mount straps (OPTION)
- I - Bed frame attachment (OPTION)

III HOW TO USE THE HANDCONTROL

Description

KEYBOARD

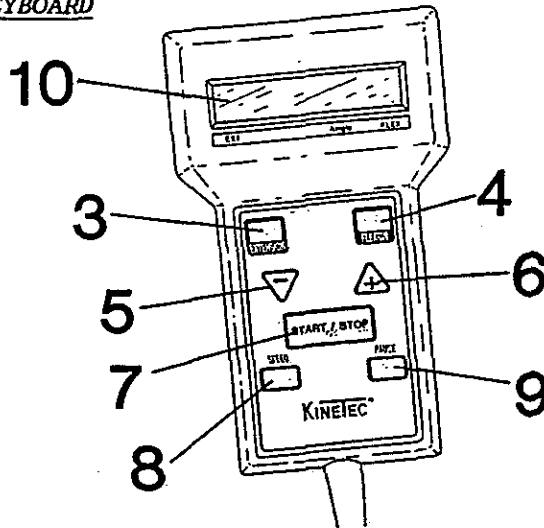- 3- EXTENSION key
- 4- FLEXION key
- 5- Decrease key
- 6- Increase key
- 7- START / STOP key
- 8- SPEED key
- 9- PAUSE key
- 10- DISPLAY

Fig.3

DISPLAY

The liquid crystal display shows up to 16 characters.


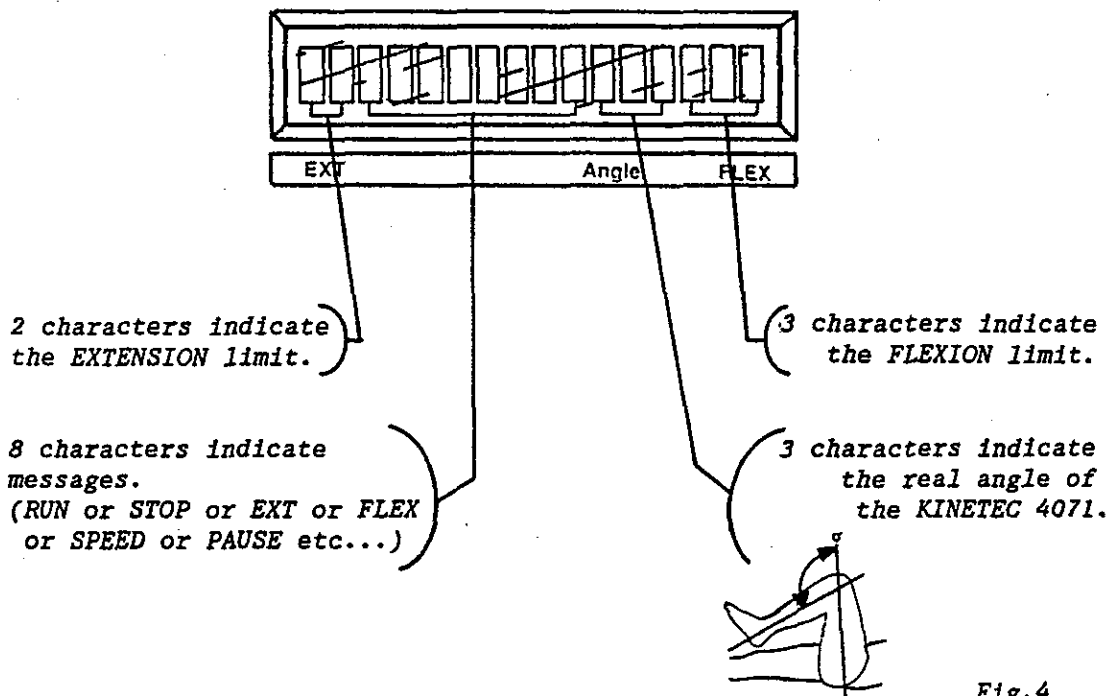
Fig.4

Adjusting the EXTENSION limit

Check that the switch (11) is in position

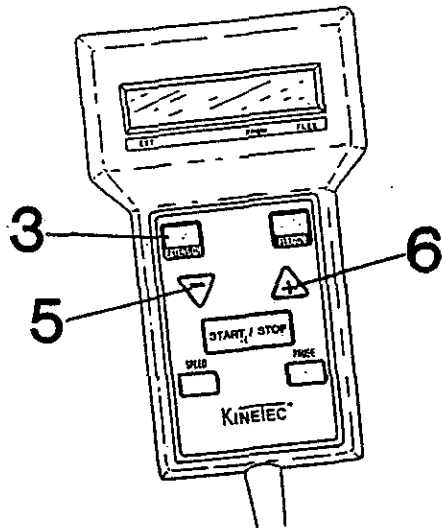Modification of the **EXTENSION** limit can be done in either the running or stopped position.

- Press the key

The display indicates

- Press  or 

until the display indicates the desired **EXTENSION** limit. *

- MINI : -9°

- MAXI : 99°


Adjusted in 1° increments.

NOTE

After 3 seconds the KINETEC 4071 will beep and the display will flash once to confirm that one of the operation settings has been modified.

Adjusting the FLEXION limit

Check that the switch (11) is in position

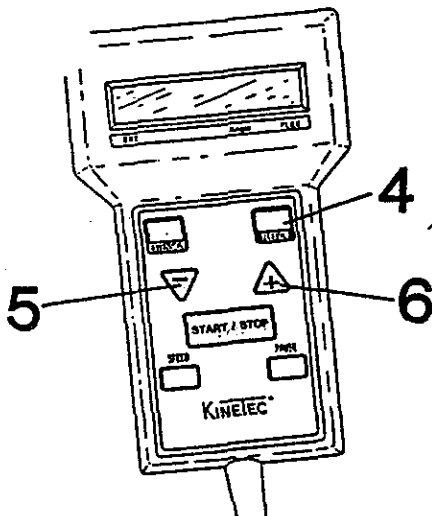
Modification of the **FLEXION** limit can be done in either the running or stopped position.

- Press the key

The display indicates

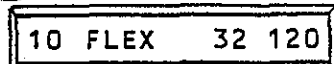- Press  or 

until the display indicates the desired **FLEXION** limit. *

- MINI : -4°

- MAXI : 120°

Adjusted in 1° increments.

NOTE


After 3 seconds the KINETEC 4071 will beep and the display will flash once to confirm that one of the operation settings has been modified.

* The legend "5° LIMIT" will appear to indicate that the flexion limit is always 5° greater than the extension limit.

MANUAL use

This function allows the limits of movement tolerated by the patient to be measured, and possibly stored.

Check that the switch (11) is in position

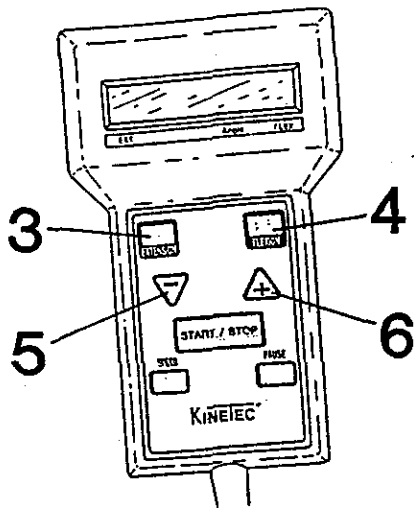Stop the KINETEC 4071

Hold the key for at least 3 seconds :

- on the  key for EXTENSION movement

- on the  key for FLEXION movement

The display indicates

10 MANUAL 15 100

On releasing the key the movement stops.

If there is no further use of the

or  keys for 3 seconds the manual mode switches off

NOTE

10 STOP 20 100

The use of the manual function does not change the movement limits which had already been selected. To record one of the limits, press within the 3 seconds,


press the key for the direction concerned,  EXTENSION or  FLEXION.

The speed in manual mode is constant at 70°/minute.

BY-PASSING use

This function allows the pre-selected FLEXION/EXTENSION limits to be over-ridden, which allows the KINETEC 4071 to be adjusted to the maximum ranges tolerated by the patient, without having stopped the machine, but with complete safety.

Check the switch (11) is in position

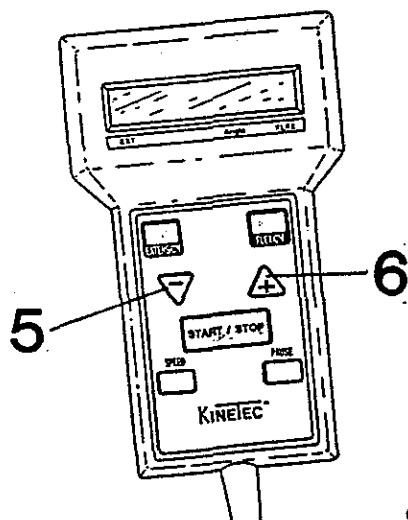Start the KINETEC 4071

Maintain the hold for 3 seconds :

- on the  key to exceed the EXTENSION limit seen on the display.

- on the  switch to exceed the FLEXION limit seen on the display.

The display shows

10 ByPass 2 90

On releasing the switch the movement reverses.

When the angle of the patient's leg exceeds the pre-selected EXTENSION or FLEXION limit, the displays flashes

7 ByPass 7 90

or

10 ByPass 101 101

Release the hold on the  or  keys when the patient doesn't tolerate greater FLEXION or EXTENSION and the movement reverses.

To record the limits defined by the BY-PASS mode, press within the 3 seconds on the switch for the direction concerned. If not, after 3 seconds the BY-PASS mode switches off without changing the limits of movements.

IV USING THE KINETEC 4071

Set up of the pads on the KINETEC 4071

The pads are numbered

- 4071.1 : for the thigh
- 4071.2 : for the leg
- 4071.3 : for the foot
- 4071.4 : for the ankle

Number for order

the complet pads : 50000183.

Set up the pad 4071.1

Make sure the part (P1) is put under the knob (B1)(Fig.A).

Set up the foot pad 4071.3

Set up the foot pad on the sole (E).

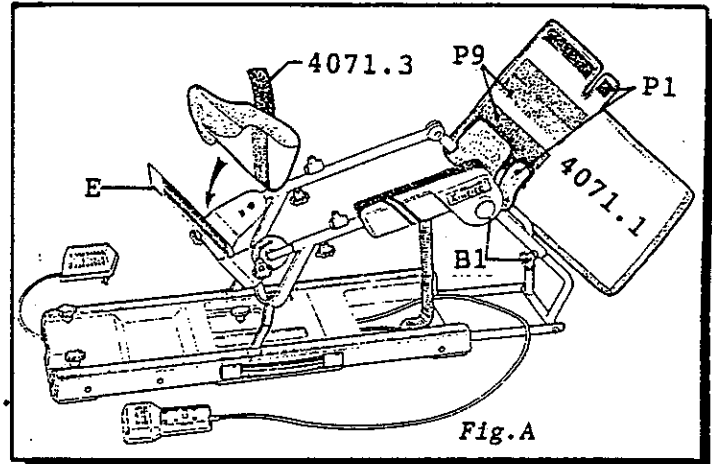
Fig.A

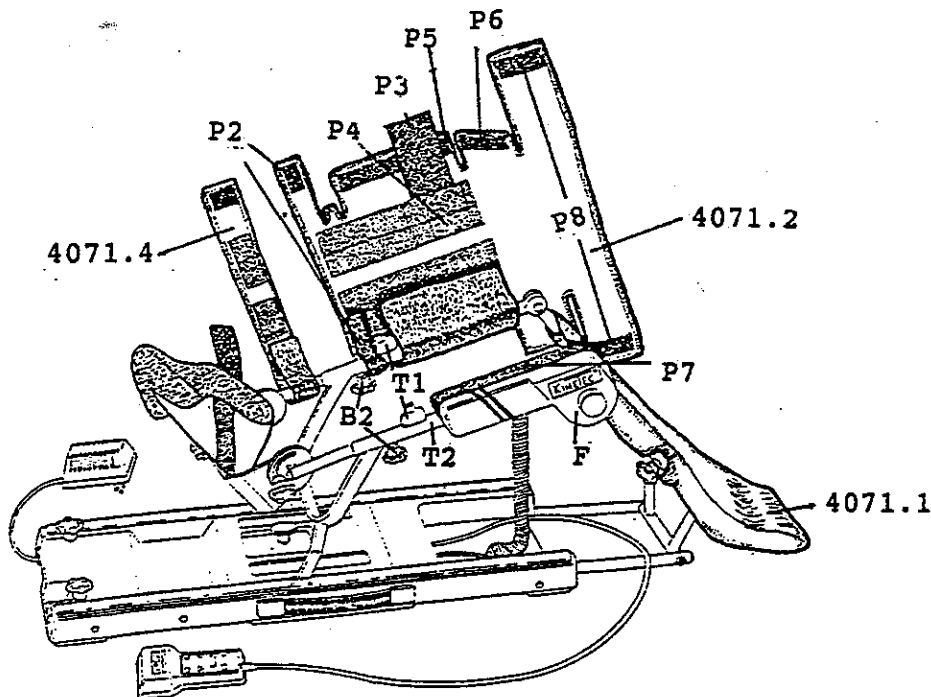
Fig.B

Set up the pad 4071.2 (Fig.B)

- The pad 4071.2 covers partially the pad 4071.1.
- Make sure the 2 parts (P2) are putting between the 2 knobs (B2) and the 2 tubes (T1).
- The part (P3) is glide between the tube (T2) of the leg support and the motor block (F), it is fixing on (P4).
- The part (P5) and (P6) are fixing on the motor cover on (P7). Make sure the part (P6) is covering the motor cover angle (F).
- The 2 parts (P8) are fixing under the thigh pad 4071.1 on (P9)(Fig.A).


Set up the ankle pad 4071.4 (Fig.B)

Using only the ankle pad 4071.4 for a long leg.
It covers the foot support.

TRANSPORT position

The KINETEC 4071 is "TELESCOPIC", which allows it to be shortened for transport.

- Putting the KINETEC 4071 in TRANSPORT position.

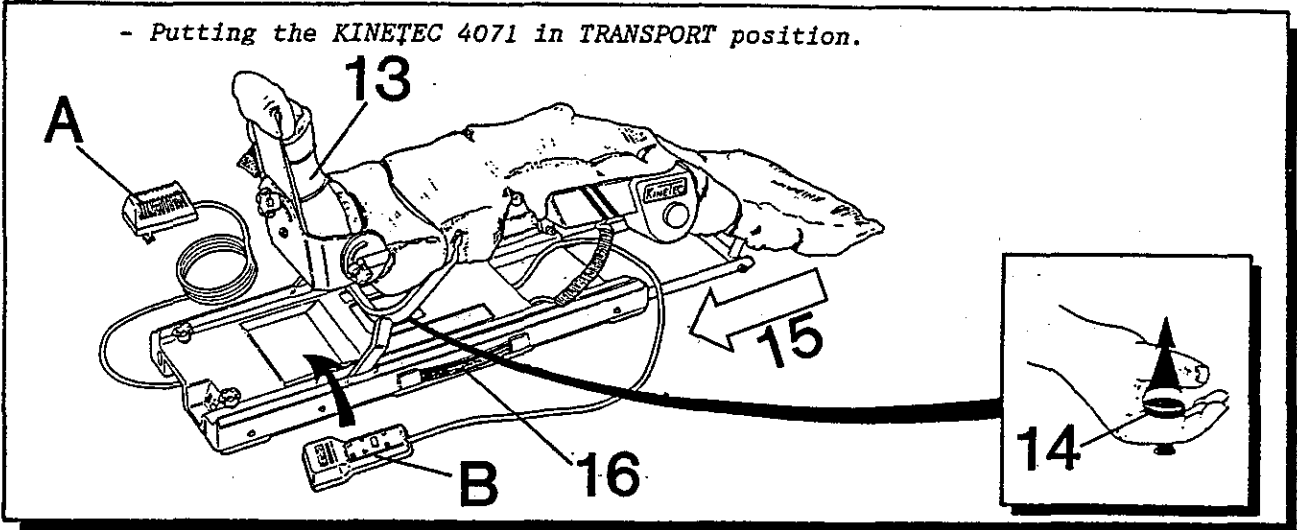
- Adjust the KINETEC 4071 to its smallest size.(mini femur)
- Stop the machine at 18° of flexion.
- Switch off and unplug the transformer (A).
- Put the transformer (A) within the footpad (13) and wrap the cable around the footplate.
- Fold the footplate flat.
- Pull button (14).
- Slide the leg support (15) in the direction of arrow.
- Attach the handcontrol (B).(Velcro)
- Always use the carrying handle (16) for lifting the KINETEC 4071.

IMPORTANT

THE KINETEC 4071 WILL NOT RUN IN THE TRANSPORT POSITION.

WORKING position


- Putting the KINETEC 4071 in the WORKING position.

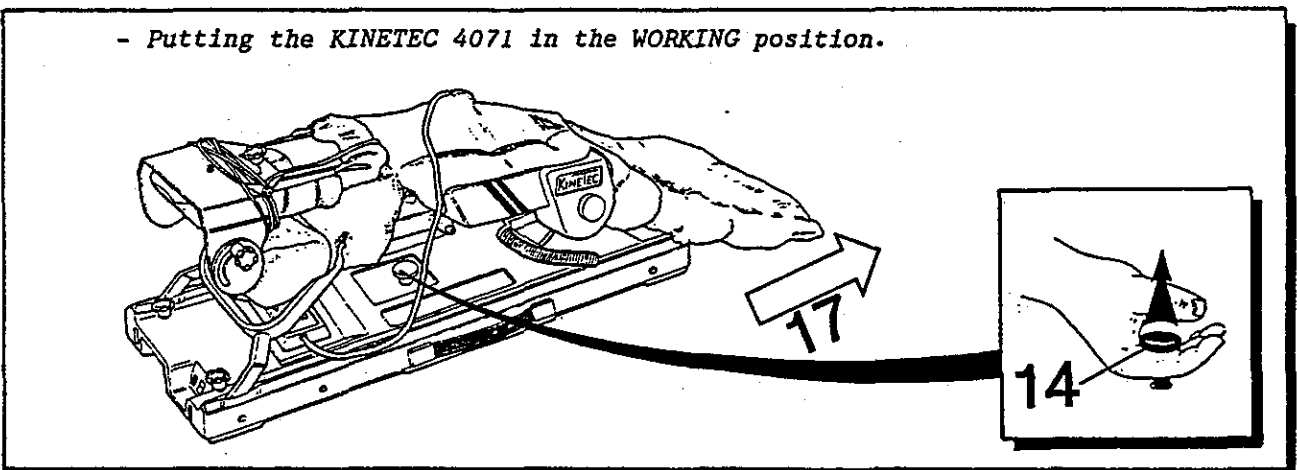- Pull button (14).
- Slide the leg supports (17) out in the direction of the arrow.


Positionning of the foot plate

The articulated footplate allows the positioning of the foot from 40° in PLANTAR FLEXION to 20° in DORSIFLEXION and 30° in EXTERNAL or INTERNAL ROTATION.

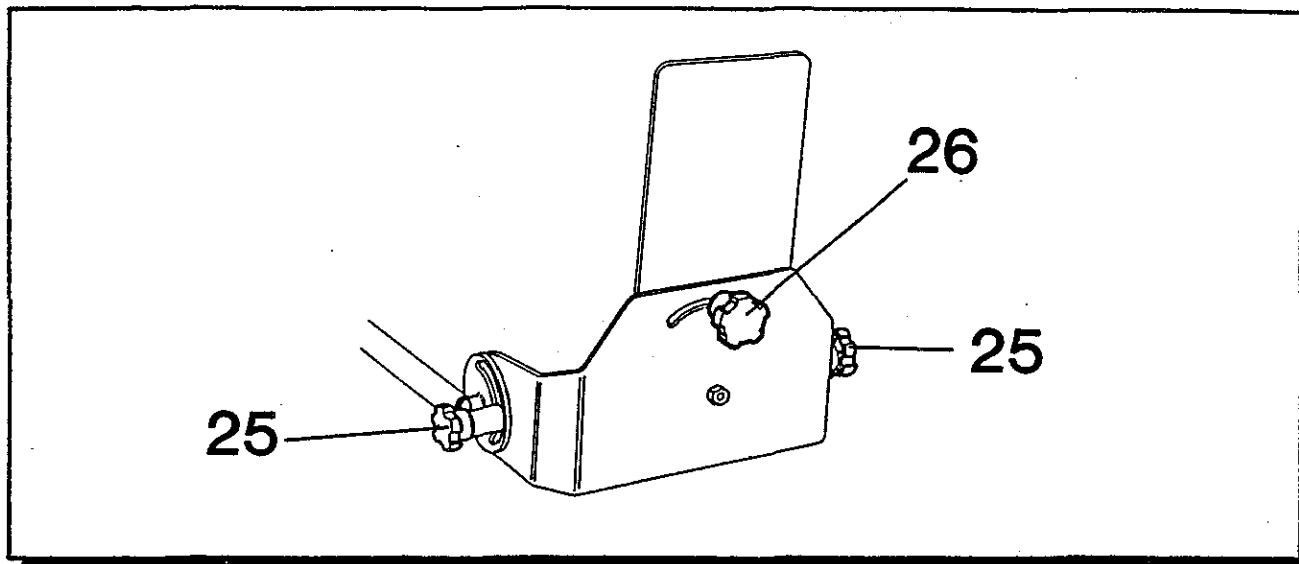1. To get the foot in DORSIFLEXION or PLANTAR FLEXION :
 - Loosen the 2 knobs (25).
 - Set the position of the foot.
 - Tighten the 2 knobs (25).
2. To get the foot in EXTERNAL or in INTERNAL ROTATION :
 - Loosen the knob (26).
 - Set the position of the foot.
 - Tighten the knob (26).

Protection of the opposite foot

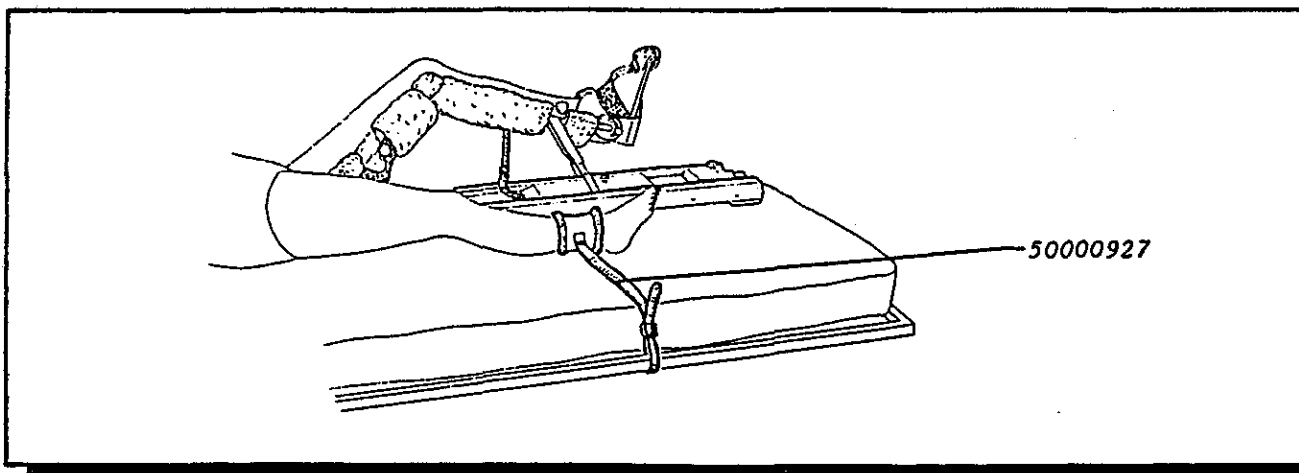EXTRA SAFETY

When using the machine in bed it is recommended to restrain the opposite foot so that it does not impede the correct working of the KINETEC 4071. Use the strap supplied with the machine for this purpose. Order number for a replacement strap : 50000927.

VI INSTRUCTIONS IN CASE OF A BREAKDOWN

After switching on the power supply to the KINETEC 4071 using the switch POWER (2)(cf. Page 3) :

- The KINETEC 4071 fails to operate, and if the handcontrol does not give any information :
 - Check the working position.
 - Check the electrical connections.
 - Check at the mains socket that power supply is operational by using another electrical device.
- The KINETEC 4071 fails to operate, but the display indicate

15 STOP	25 120
---------	--------

 - Press the START / STOP key once again.
- The KINETEC 4071 still fails to operate after that :
 - Return the machine to your nearest specialist.

GUARANTEE LIMITS

The KINETEC guarantee is strictly limited to free replacement or to repairs at the plant of the part or parts acknowledged as defective.

KINETEC guarantees its passive motion device for 1 year from the date of purchase by the consumer against any manufacturing defect.

KINETEC alone is qualified to judge as to the application of the guarantee on its device.

The guarantee cannot be applied if the apparatus has been subject to abnormal use or has been used in circumstances other than those stipulated in our user manual.

The guarantee will not be applied, either, in case of deterioration or accident due to negligence, lack of supervision or maintenance resulting from a transformation of the equipment or from an attempt to repair the equipment.

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